

## **CheckLite™ HS Set**

Luciferin-luciferase ATP assay reagent kit for hypersensitive viable count assay.



## Characteristics of kit

- CheckLite<sup>™</sup> HS Set is a kit for microbial biomass assays based on the ATP measuring method.
- CheckLite™ contains thermostable firefly luciferase and ATP eliminating enzymes.
- CheckLite<sup>™</sup> is characterized by its ease of use, rapid performance and high sensitivity.
- CheckLite™ can detect 1,000 cells of coliforms/ml (100 cell/assay).

## Composition of kit

Luciferin-luciferase reagent HS (lyophilized)
Reconstitution buffer (5.5 ml)
ATP eliminating reagent (lyophilized)
Reconstitution buffer (5.5 ml)
ATP releasing reagent (5.5 ml)


## (100 assays/kit)

2

2

2

2

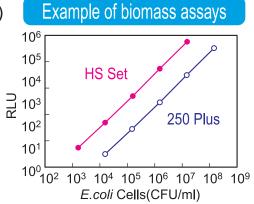

CheckLite<sup>™</sup> HS Set should be stored at 2-8°C in a refrigerator. DO NOT FREEZE.

Manufacturer, Contact



Kikkoman Biochemifa Company

URL: http://biochemifa.kikkoman.co.jp/

2-1-1, Nishi-shinbashi, Minato-ku, Tokyo 105-0003, Japan Phone: +81-3-5521-5492 Fax: +81-3-5521-5498 E-mail: biochemifa@mail.kikkoman.co.jp